



**Expanded Range
High Sensitivity**

SALIVARY CORTISOL

ENZYME IMMUNOASSAY KIT

For Research Use Only
Not for use in Diagnostic Procedures

Item No. 1-3002, (Single) 96-Well Kit;
1-3002-5, (5-Pack) 480 Wells



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

TABLE OF CONTENTS

Intended Use	3
Introduction	3
Test Principle	4
Safety Precautions	4
General Kit Use Advice	5
Storage	5
pH Indicator	5
Specimen Collection	6
Sample Handling and Preparation	6
Materials Supplied with Single Kit	7
Materials Needed But Not Supplied	8
Reagent Preparation	8
Procedure	9
Quality Control	10
Calculations	10
Typical Results	11
Limitations	12
Salivary Cortisol Example Ranges	13
Salivary Cortisol EIA Kit Performance Characteristics	14
References	19
Seller's Limited Warranty	20



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Intended Use

The Salimetrics® Cortisol Enzyme Immunoassay Kit is a competitive immunoassay specifically designed and validated for the quantitative measurement of salivary Cortisol. It is not intended for diagnostic use. It is intended only for research use in humans and some animals. Salimetrics has not validated this kit for serum or plasma samples.


Please read the complete kit insert before performing this assay. Failure to follow kit procedure and recommendations for saliva collection and sample handling may result in unreliable values.

For further information about this kit, its application, or the procedures in this insert, please contact the technical service team at Salimetrics or your local sales representative.

Introduction

Cortisol (hydrocortisone, Compound F) is the major glucocorticoid produced in the adrenal cortex. Cortisol production has a circadian rhythm (3). Levels peak in the early morning and drop to the lowest concentration at night (4). Levels rise independently of circadian rhythm in response to stress (5). Increased Cortisol production is associated with Cushing's syndrome, while decreased Cortisol production is associated with adrenal insufficiency (e.g., Addison's disease) (6).

In blood, only about 5-10% of Cortisol is in its unbound or biologically active form. The remaining Cortisol is bound to serum proteins (7). Unbound serum Cortisol enters the saliva via intracellular mechanisms; in saliva, the majority of Cortisol remains unbound to protein (8). Salivary Cortisol levels are unaffected by salivary flow rate and are relatively resistant to degradation from enzymes or freeze-thaw cycles (8,9). Studies consistently report high correlations between serum and saliva Cortisol, indicating that salivary Cortisol levels reliably estimate serum Cortisol levels (10-12).

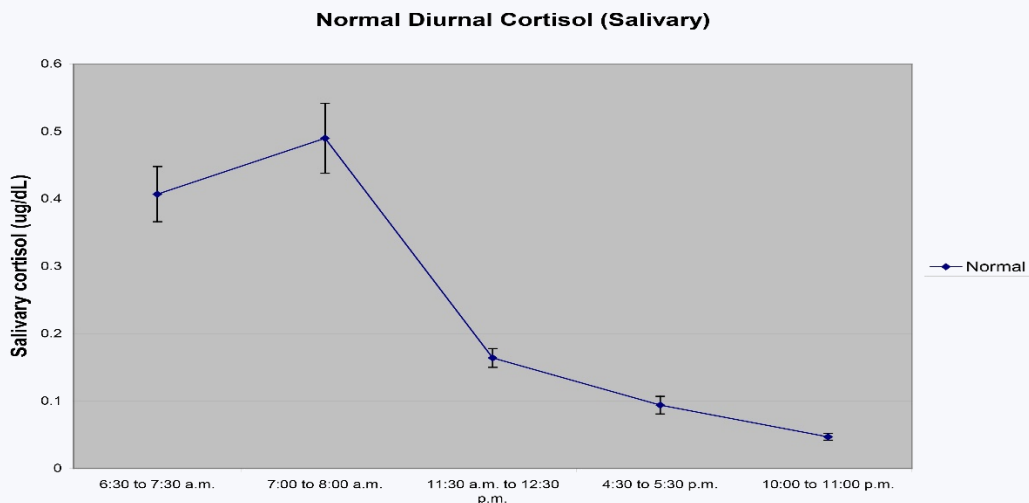

(Internal Salimetrics Data, n=26. Time of Cortisol peak will vary in individuals relative to their normal wake-up time.)



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Test Principle

This is a competitive immunoassay kit. Cortisol in standards and samples compete with Cortisol conjugated to horseradish peroxidase for the antibody binding sites on a microtitre plate. After incubation, unbound components are washed away. Bound Cortisol Enzyme Conjugate is measured by the reaction of the horseradish peroxidase enzyme to the substrate tetramethylbenzidine (TMB). This reaction produces a blue color. A yellow color is formed after stopping the reaction with an acidic solution. The optical density is read on a standard plate reader at 450 nm. The amount of Cortisol Enzyme Conjugate detected is inversely proportional to the amount of Cortisol present in the sample (13).

Safety Precautions

Read Safety Data Sheets before handling reagents.

Hazardous Ingredients

Liquid Stop Solution is caustic; use with care. We recommend the procedures listed below for all kit reagents.

Handling

Follow good laboratory practices when handling kit reagents. Laboratory coats, gloves, and safety goggles are recommended. Wipe up spills using appropriate absorbent materials while wearing protective clothing. Follow local regulations for disposal.

Emergency Exposure Measures

In case of contact, immediately wash skin or flush eyes with water for 15 minutes. Remove contaminated clothing. If inhaled, remove individual to fresh air. If individual experiences difficulty breathing call a physician.

The above information is believed to be accurate but is not all-inclusive. This information should be used only as a guide. Salimetrics will not be liable for accidents or damage resulting from misuse of product.

Safety Data Sheets are available by contacting Salimetrics at support@salimetrics.com (See www.salimetrics.com for alternative contact options).



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

General Kit Use Advice

- This kit uses break-apart microtitre strips. You may run less than a full plate. Unused wells must be stored at 2-8°C in the foil pouch with desiccant and used in the frame provided.
- Avoid microbial contamination of opened reagents. Salimetrics recommends using opened reagents within one month. Store all reagents at 2-8°C.
- The quantity of reagent provided with a single kit is sufficient for three partial runs. The volumes of wash buffer and enzyme conjugate prepared for assays using less than a full plate should be scaled down accordingly, keeping the same dilution ratio.
- Do not mix components from different lots of kits.
- To ensure highest quality assay results, pipetting of samples and reagents must be done as quickly as possible (without interruption) across the plate. Ideally, the process should be completed within 20 minutes or less.
- When using a multichannel pipette to add reagents, always follow the same sequence when adding all reagents so that the incubation time is the same for all wells.
- When running multiple plates, or multiple sets of strips, a standard curve must be run with each individual plate and/or set of strips.
- The temperature of the laboratory may affect assays. Salimetrics' kits have been validated at 68-74°F (20-23.3°C). Higher or lower temperatures may affect OD values.
- Routine calibration of pipettes and other equipment is critical for the best possible assay performance.
- When mixing plates during assay procedures, avoid speeds that spill the contents of the wells.
- We recommend saving all reagents until data analysis has confirmed a successful run to facilitate troubleshooting if necessary.
- Prior to sample addition, please label each strip to assure plate orientation and sample order when data is acquired on plate reader.

Storage

All unopened components of this kit are stable at 2-8°C until the kit's expiration date.

pH Indicator

Cortisol values from samples with a pH ≤ 3.5 or ≥ 9.0 may be inaccurate. A pH indicator in the Assay Diluent alerts the user to samples with high or low pH values. Upon addition of the Assay Diluent, acidic samples will turn yellow and alkaline samples will turn purple. Dark yellow or purple wells indicate that a pH value for that sample should be obtained using pH strips. Samples with a pH ≤ 3.5 or ≥ 9.0 should be recollected (14).



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Specimen Collection

Avoid sample collection within 60 minutes after eating a major meal or within 12 hours after consuming alcohol. Acidic or high sugar foods can compromise assay performance by lowering sample pH and influencing bacterial growth. To minimize these factors, rinse mouth thoroughly with water 10 minutes before sample is collected.

Collect whole saliva by unstimulated passive drool. Donors may tilt the head forward, allowing the saliva to pool on the floor of the mouth, then pass the saliva through the SalivaBio Collection Aid (SCA) into a polypropylene vial. Collection protocols/methods are available online at www.salimetrics.com or upon request.

Samples visibly contaminated with blood should be recollected. Samples may be screened for possible blood contamination (15,16) using our Blood Contamination EIA Kit (Item Nos. 1-1302/1-1302-5). Do not use dipsticks, which result in false positive values due to salivary enzymes.

It is important to record the time and date of specimen collection when samples are obtained due to the diurnal variation in Cortisol levels. Samples for Cushing's diagnosis should be collected at 11:00 pm.

Sample Handling and Preparation

After collection it is important to keep samples cold in order to avoid bacterial growth in the specimen. Refrigerate sample within 30 minutes, and freeze at or below -20°C within 4 hours of collection. (Samples may be stored at -20°C for up to 6 months.) For long term storage, refer to the Salimetrics Collection and Handling Advice Booklet.

Do not add sodium azide to saliva samples as a preservative, as it may cause interference in the immunoassay.

On day of assay, thaw the saliva samples completely, vortex, and centrifuge at 1500 x g for 15 minutes. Freezing saliva samples will precipitate mucins. Centrifuging removes mucins and other particulate matter which may interfere with antibody binding and affect results. Samples should be at room temperature before adding to assay plate. Pipette clear sample into appropriate wells. Re-freeze saliva samples as soon as possible after adding to the assay plate. Re-centrifuge saliva samples each time that they are thawed. Avoid multiple freeze-thaw cycles.



Materials Supplied with Single Kit

	Item	Quantity/Size
1	Microtitre Plate Coated with monoclonal anti-Cortisol antibodies.	1/96 well
2	Cortisol Standard In a saliva-like matrix. Ready to use, traceable to NIST standard: 3.0, 1.0, 0.333, 0.111, 0.037, 0.012 µg/dL (82.77, 27.59, 9.19, 3.06, 1.02, 0.33 nmol/L). Contains: Cortisol, buffer, preservative.	6 vials / 500 µL each
3	Cortisol Controls High, Low, in a saliva-like matrix. Ready to use. Contain: Cortisol, buffer, preservative.	2 vials / 500 µL each
4	Cortisol Enzyme Conjugate Concentrate. Dilute before use with Assay Diluent. (See step 5 of Procedure.) Contains: Cortisol conjugated to HRP, preservative.	1 vial / 50 µL
5	Assay Diluent Contains: phosphate buffer, pH indicator, preservative.	1 bottle / 60 mL
6	Wash Buffer Concentrate (10X) Dilute before use according to Reagent Preparation. Contains: phosphate buffer, detergent, preservative.	1 bottle / 100 mL
7	TMB Substrate Solution Non-toxic, ready to use.	1 bottle / 25 mL
8	Stop Solution	1 bottle / 12.5 mL



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Materials Needed But Not Supplied

- Precision pipette to deliver 15 and 25 μL
- Precision multichannel pipette to deliver 50 μL and 200 μL
- Vortex
- Plate rotator with 0.08-0.17 inch orbit capable of 500 rpm
- Plate reader with 450 nm and 490 to 492 nm reference filters
- Computer software for data reduction
- Deionized water
- Reagent reservoirs
- One disposable polypropylene tube to hold at least 24 mL
- Pipette tips
- Serological pipette to deliver up to 24 mL
- Centrifuge capable of 1500 x g

Reagent Preparation

- Bring all reagents to room temperature and mix before use. A minimum of 1.5 hours is recommended for the 24 mL of Assay Diluent used in Step 5 (conjugate dilution) to come to room temperature.
- Bring Microtitre Plate to room temperature before use. ***It is important to keep the foil pouch with the plate strips closed until warmed to room temperature, as humidity may have an effect on the coated wells.***
- Prepare 1X wash buffer by diluting Wash Buffer Concentrate (10X) 10-fold with room-temperature deionized water (100 mL of Wash Buffer Concentrate (10X) to 900 mL of deionized H_2O). ***Dilute only enough for current day's use and discard any leftover reagent.*** (If precipitate has formed in the concentrated wash buffer, it may be heated to 40°C for 15 minutes. Cool to room temperature before use in assay.)



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Procedure

Step 1: Read and prepare reagents according to the Reagent Preparation section before beginning assay. Determine your plate layout. Here is a suggested layout. (Standards, controls, and saliva samples should be assayed in duplicate.)

	1	2	3	4	5	6	7	8	9	10	11	12
A	3.000 Std	3.000 Std	Ctrl-L	Ctrl-L								
B	1.000 Std	1.000 Std	SMP-1	SMP-1								
C	0.333 Std	0.333 Std	SMP-2	SMP-2								
D	0.111 Std	0.111 Std	SMP-3	SMP-3								
E	0.037 Std	0.037 Std	SMP-4	SMP-4								
F	0.012 Std	0.012 Std	SMP-5	SMP-5								
G	0 Std	0 Std	SMP-6	SMP-6								
H	Ctrl-H	Ctrl-H	SMP-7	SMP-7								

Step 2: Keep the desired number of strips in the strip holder and place the remaining strips back in the foil pouch with the desiccant. Store sealed at 2-8°C.

Caution: Do not insert wells from one plate into a different plate

Step 3: Pipette 24 mL of Assay Diluent into the disposable tube. (Scale down proportionally if using less than the entire plate.) Set aside for Step 5.

Step 4:

- Pipette 25 µL of standards, controls, and saliva samples into appropriate wells.
- Pipette 25 µL of Assay Diluent into 2 wells to serve as the zero.

Step 5: Dilute the Enzyme Conjugate 1:1600 by adding 15 µL of the conjugate to the 24 mL tube of Assay Diluent. (Scale down proportionally if not using the entire plate.) Conjugate tube may be centrifuged for a few minutes to bring the liquid down to the tube bottom. Immediately mix the diluted conjugate solution and add 200 µL to each well using a multichannel pipette.

Step 6: Mix plate on a plate rotator for 5 minutes at 500 rpm and incubate at room temperature for a total of 1 hour.



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Step 7: Wash the plate 4 times with 1X wash buffer. A plate washer is recommended. However, washing may be done by gently squirting wash buffer into each well with a squirt bottle, or by pipetting 300 μL of wash buffer into each well and then discarding the liquid over a sink. After each wash, the plate should be thoroughly blotted on paper towels before turning upright. If using a plate washer, blotting is still recommended after the last wash.

Step 8: Add 200 μL of TMB Substrate Solution to each well with a multichannel pipette.

Step 9: Incubate the plate in the dark (covered) at room temperature for 30 minutes, mixing for 5 minutes on a plate rotator at 500 rpm.

Step 10: Add 50 μL of Stop Solution with a multichannel pipette.

Step 11:

- Mix on a plate rotator for 3 minutes at 500 rpm. If green color remains, continue mixing until green color turns to yellow. Be sure all wells have turned yellow.

Caution: Spillage may occur if mixing speed exceeds 600 rpm.

- Wipe off bottom of plate with a water-moistened, lint-free cloth and wipe dry.
- Read in a plate reader at 450 nm. Read plate within 10 minutes of adding Stop Solution. (For best results, a secondary filter correction at 490 to 492 nm is recommended.)

Quality Control

The Salimetrics' High and Low Cortisol Controls should be run with each assay. The control ranges established at Salimetrics are to be used as a guide. Each laboratory should establish its own range. Variations between laboratories may be caused by differences in techniques and instrumentation.

Calculations

1. Compute the average optical density (OD) for all duplicate wells.
2. Calculate the percent bound (B/B_0) for each standard, control, and saliva sample by dividing the OD of each well (B) by the average OD for the zero (B_0).
3. Determine the concentrations of the controls and saliva samples by interpolation using data reduction software. We recommend using a 4-parameter non-linear regression curve fit.
4. Samples with Cortisol values greater than 3.0 $\mu\text{g/dL}$ (82.77 nmol/L) should be diluted with Assay Diluent and rerun for accurate results. If a dilution of the sample is used, multiply the assay results by the dilution factor.

A new Standard Curve must be run with each full or partial plate.




101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Typical Results

The results shown below are for illustration only and should not be used to calculate results from another assay.

Well	Standard	Average OD (B)	B/B ₀	Cortisol (µg/dL)
A1,A2	S1	0.094	0.079	3.000
B1,B2	S2	0.236	0.194	1.000
C1,C2	S3	0.524	0.417	0.333
D1,D2	S4	0.897	0.671	0.111
E1,E2	S5	1.219	0.847	0.037
F1,F2	S6	1.379	0.941	0.012
G1,G2	Zero	1.496	1	0

Example: HS Cortisol 4-Parameter Curve Fit

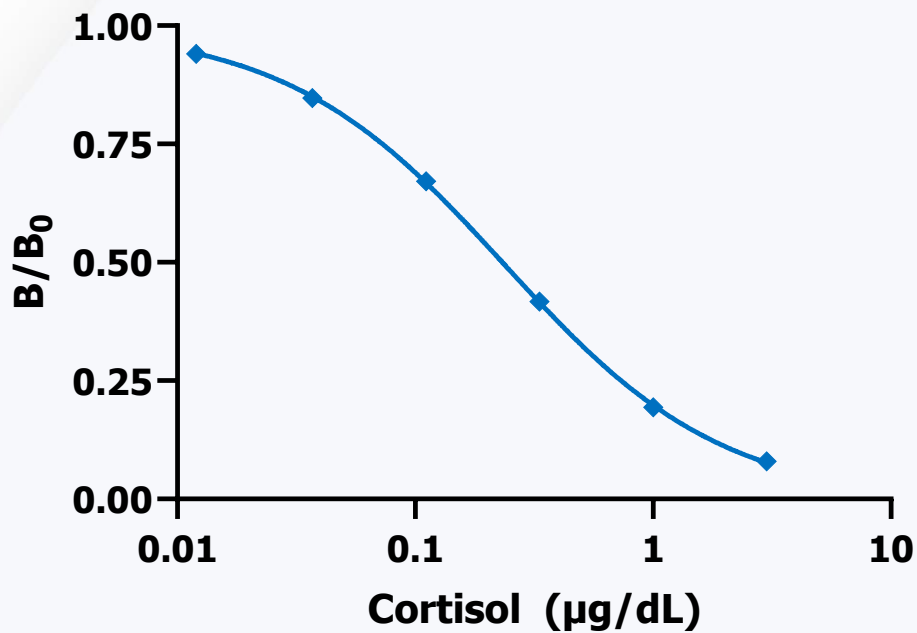

101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Limitations

- Diagnosis of Cushing's syndrome should be confirmed by additional diagnostic tests for the disease, such as low-dose dexamethasone suppression testing.
- Cortisol levels are elevated during the later stages of pregnancy and in women on oral contraceptives or after long-term use of oral contraceptives (12,17).
- Some studies show developmental differences in Cortisol as well as an association between Cortisol and weight (18).
- Elevated Cortisol levels can be found in conditions of sepsis, infection, chronic liver disease, and renal failure. Low Cortisol levels result from liver disease, pituitary hyposecretion, hypothyroidism, or steroid therapy.
- Samples with Cortisol values greater than 3.0 µg/dL (82.77 nmol/L) should be diluted with Assay Diluent and rerun for accurate results. To obtain the final Cortisol concentration, multiply the concentration of the diluted sample by the dilution factor.
- A pH value should be obtained on samples that appear yellow or purple after the diluted conjugate solution is added and the plate is mixed (Step 6). Samples with pH values ≤ 3.5 or ≥ 9.0 should be recollected.
- See "Specimen Collection" recommendations to ensure proper collection of saliva specimens and to avoid interfering substances.
- Samples collected with sodium azide are unsuitable for this assay.
- Any quantitative results indicating abnormal Cortisol levels should be followed by additional testing and evaluation.



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Salivary Cortisol Example Ranges*

Group	Number	Overall Range (µg/dL)
Children, neonatal	275	ND - 3.417
Children, age 6 months	165	ND - 2.734

Group	Number	23:00 hrs (µg/dL)
Normal subjects	19	0.007 – 0.115
Cushing's subjects	21	0.130 – 2.972

Group	Number	AM Range (µg/dL)	PM Range (µg/dL)
Children, ages 2.5-5.5	112	0.034 - 0.645	0.053 - 0.607
Children, ages 8-11	285	0.084 - 0.839	ND - 0.215
Adolescents, ages 12-18	403	0.021 - 0.883	ND - 0.259
Adult males, ages 21-30	26	0.112 - 0.743	ND - 0.308
Adult females, ages 21-30	20	0.272 - 1.348	ND - 0.359
Adult males, ages 31-50	67	0.122 - 1.551	ND - 0.359
Adult females, ages 31-50	31	0.094 - 1.515	ND - 0.181
Adult males, ages 51-70	28	0.112 - 0.812	ND - 0.228
Adult females, ages 51-70	23	0.149 - 0.739	0.022 - 0.254
All adults	192	0.094 - 1.551	ND - 0.359

*To be used as a guide only. Each laboratory should establish its own range.

ND = None detected

Expected ranges for neonates to 5.5 years were derived using the Salimetrics Salivary Cortisol Immunoassay Kit.

Expected ranges for 8 to 18 years were reported from an unpublished manuscript, Pennsylvania State University's Behavioral Endocrinology Laboratory. Adult ranges were obtained from published literature (7).




101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

HS Salivary Cortisol EIA Kit Performance Characteristics

Linearity of Assay

Saliva Sample	Dilution	Expected (µg/dL)	Observed (µg/dL)	Recovery (%)
1		0.07	0.07	100%
	1:9	0.36	0.34	108%
	2:8	0.63	0.61	104%
	3:7	0.93	0.88	106%
	4:6	1.13	1.15	98%
	5:5	1.45	1.42	102%
	6:4	1.64	1.69	97%
	7:3	1.88	1.96	96%
	8:2	2.27	2.23	102%
	9:1	2.49	2.50	99%
2		2.77	2.77	100%

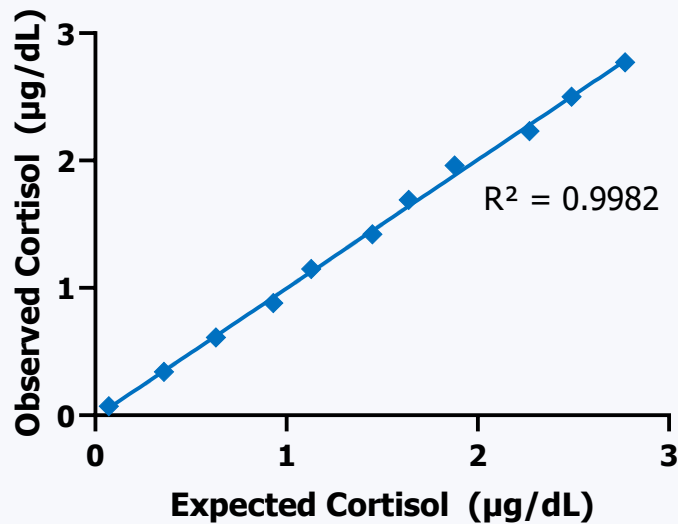

101 Innovation Boulevard • Suite 302 • State College, PA 16803
 1.800.790.2258 • support@salimetrics.com • salimetrics.com

Spike and Recovery

Five saliva samples containing different levels of an endogenous Cortisol were spiked with known quantities of Cortisol and assayed.

Saliva Sample	Endogenous (µg/dL)	Added (µg/dL)	Expected (µg/dL)	Observed (µg/dL)	Recovery (%)
1	0.071	2.00	2.07	2.20	106%
2	0.071	0.20	0.27	0.28	104%
3	0.071	0.04	0.11	0.11	98%
4	0.078	2.33	2.41	2.33	97%
5	0.078	0.20	0.28	0.31	113%
6	0.080	0.04	0.12	0.12	103%
7	0.860	0.20	1.06	1.16	109%
8	0.890	0.04	0.93	1.02	109%



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Sample Dilution Recovery

Four saliva samples were diluted with Assay Diluent and assayed.

Saliva Sample	Dilution Factor	Expected (µg/dL)	Observed (µg/dL)	Recovery (%)
1	undiluted	N/A	0.73	N/A
	1:2	0.37	0.39	107%
	1:4	0.18	0.20	111%
	1:8	0.09	0.10	111%
	1:16	0.05	0.05	105%
2	undiluted	N/A	0.80	N/A
	1:2	0.40	0.40	101%
	1:4	0.20	0.19	97%
	1:8	0.10	0.09	94%
	1:16	0.05	0.05	110%
3	undiluted	N/A	0.61	N/A
	1:2	0.31	0.30	98%
	1:4	0.15	0.15	101%
	1:8	0.08	0.08	108%
	1:16	0.04	0.04	108%
4	undiluted	N/A	2.89	N/A
	1:2	1.45	1.53	105%
	1:4	0.72	0.77	106%
	1:8	0.36	0.42	115%
	1:16	0.18	0.20	108%

Precision

The intra-assay precision was determined from the mean of 20 replicates each.

Saliva Sample	N	Mean (µg/dL)	Standard Deviation (µg/dL)	Coefficient of Variation (%)
1	20	2.07	0.08	4%
2	20	1.14	0.05	4%
3	20	0.42	0.01	3%
4	20	0.16	0.01	5%
5	20	0.06	0.00	7%

The inter-assay precision was determined from the mean of average duplicates for 20 separate runs.

Saliva Sample	N	Mean (µg/dL)	Standard Deviation (µg/dL)	Coefficient of Variation (%)
1	20	1.99	0.05	3%
2	20	1.16	0.05	4%
3	20	0.43	0.01	3%
4	20	0.18	0.01	9%
5	20	0.06	0.01	11%

Sensitivity

Analytical Sensitivity

The lower limit of sensitivity was determined by interpolating the mean optical density minus 2 SDs of 10 sets of duplicates at the 0 µg/dL level. The minimal concentration of Cortisol that can be distinguished from 0 is **0.007 µg/dL**.

Functional Sensitivity

The functional sensitivity was determined by assaying 60 samples at a concentration level resulting in a CV of approximately 20%. The functional sensitivity of the salivary Cortisol ELISA is **0.018 µg/dL**.



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Correlation with Serum

The correlation between serum and saliva Cortisol was determined by assaying 49 matched samples using the Diagnostic Systems Laboratories Serum Cortisol EIA and the Salimetrics HS Salivary Cortisol EIA.

The correlation between saliva and serum was highly significant, $r(47) = 0.91$, $p < 0.0001$.

Antibody Specificity

Compound	Spiked Concentration (ng/mL)	% Cross-reactivity in HS Salivary Cortisol EIA
Prednisolone	100	0.568
Prednisone	1000	ND
Cortisone	1000	0.130
11-Deoxycortisol	500	0.156
21-Deoxycortisol	1000	0.041
17 α -Hydroxyprogesterone	1000	ND
Dexamethasone	1000	19.2
Triamcinolone	1000	0.086
Corticosterone	10,000	0.214
Progesterone	1000	0.015
17 β -Estradiol	10	ND
DHEA	10,000	ND
Testosterone	10,000	0.006
Transferrin	66,000	ND
Aldosterone	10,000	ND

ND = None detected (<0.004)



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

References

1. Raff, H., Raff, J.L., & Findling, J.W. (1998). Late-night salivary cortisol as a screening test for Cushing's Syndrome. *J Clin Endocrinol Metab*, *83*(8), 2681-86.
2. Papanicolaou, D.A., Mullen, N., Kyrou, I., & Nieman, L.K. (2002). Night time salivary cortisol: A useful test for the diagnosis of Cushing's Syndrome. *J Clin Endocrinol Metab*, *87*(10), 4515-21.
3. Dorn, L.D., Lucke, J.F., Loucks, T.L., & Berga, S.L. (2007). Salivary cortisol reflects serum cortisol: Analysis of circadian profiles. *Ann Clin Biochem*, *44* (Pt 3), 281-84.
4. Chernow, B., Alexander, H.R., Smallridge, R.C., Thompson, W.R., Cook, D., Beardsley, D., Fink, M.P., et al. (1987). Hormonal responses to graded surgical stress. *Arch Intern Med*, *147*(7), 1273-78.
5. Kreiger, D.T. (1975). Rhythms of ACTH and corticosteroid secretion in health and disease and their experimental modification. *J Steroid Biochem*, *6*(5), 758-91.
6. Rothfield, B. (1974). Plasma cortisol. In: B. Rothfield (ed.), *Nuclear medicine – in vitro* (pp. 120-5). Philadelphia, PA: Lippincott.
7. Aardal, E., & Holm, A. (1995). Cortisol in saliva – Reference ranges in relation to cortisol in serum. *Eur J Clin Chem Clin Biochem*, *33*(12), 927-32.
8. Vining, R.F., & McGinley, R.A. (1987). The measurement of hormones in saliva: Possibilities and pitfalls. *J Steroid Biochem*, *27*(1-3), 81-94.
9. Garde, A.H., & Hansen, A.M. (2005). Long-term stability of salivary cortisol. *Scand J Clin Lab Invest*, *65*(5), 433-36.
10. Francis, S.J., Walker, R.F., Riad-Fahmy, D., Hughes, D., Murphy, J.F., & Gray, O.P. (1987). Assessment of adrenocortical activity in term newborn infants using salivary cortisol determinations. *J Pediatr*, *111*(1), 129-33.
11. Hiramatsu, R. (1981). Direct assay of cortisol in human saliva by solid phase radio immunoassay and its clinical applications. *Clin Chim Acta*, *117*(2), 239-49.
12. Vining, R.F., McGinley, R.A., Maksvytis, J.J., & Ho, K.Y. (1983). Salivary cortisol: A better measure of adrenal cortical function than serum cortisol. *Ann Clin Biochem*, *20*(Pt 6), 329-35.
13. Chard, T. (1990). *An introduction to radioimmunoassay and related techniques* (4th ed.). Amsterdam: Elsevier.
14. Schwartz, E.B., Granger, D.A., Susman, E.J., Gunnar, M.R., & Laird, B. (1998). Assessing salivary cortisol in studies of child development. *Child Dev*, *69*(6), 1503-13.
15. Kivlighan, K.T., Granger, D.A., Schwartz, E.B., Nelson, V., & Curran, M. (2004). Quantifying blood leakage into the oral mucosa and its effects on the measurement of cortisol, dehydroepiandrosterone, and testosterone in saliva. *Horm Behav*, *46*(1), 39-46.
16. Schwartz, E., & Granger, D.A. (2004). Transferrin enzyme immunoassay for quantitative monitoring of blood contamination in saliva. *Clin Chem*, *50*(3), 654-56.
17. Susman, E.J., Schmeelk, K.H., Worrall, B.S., Granger, D.A., Ponirakis, A., & Chrousos, M.D. (1999). Corticotropin-releasing hormone and cortisol: Longitudinal associations with depression and antisocial behavior in pregnant adolescents. *J Amer Acad Child Adolesc Psychiatry*, *38*(4), 460-67.
18. Kiess, W., Meidert, A., Dressendorfer, R.A., Schriever, K., Kessler, U., Konig, A., Schwarz, H.P., & Strasburger, C.J. (1995). Salivary cortisol levels throughout childhood and adolescence: Relation with age, pubertal state, and weight. *Pediatr Res*, *37*(4 Pt 1), 502-6



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com

Seller's Limited Warranty

"Seller warrants that all goods sold hereunder will be free from defects in material and workmanship. Upon prompt notice by Buyer of any claimed defect, which notice must be sent within thirty (30) days from date such defect is first discovered and within three months from the date of shipment, Seller shall, at its option, either repair or replace the product that is proved to Seller's satisfaction to be defective. All claims should be submitted in written form. This warranty does not cover any damage due to accident, misuse, negligence, or abnormal use. Liability, in all cases, will be limited to the purchased cost of the kit.

It is expressly agreed that this limited warranty shall be in lieu of all warranties of fitness and in lieu of the warranty of merchantability. Seller shall not be liable for any incidental or consequential damages that arise out of the installation, use or operation of Seller's product or out of the breach of any express or implied warranties."

<p>Salimetrics, LLC (East Coast) 101 Innovation Blvd., Suite 302 State College, PA 16803, USA (T) 814.234.2617 (F) 814.238.1608 800-790-2258 (USA & Canada only) www.salimetrics.com support@salimetrics.com</p>	<p>Salimetrics, LLC (West Coast) 5962 La Place Court, Suite 275 Carlsbad, CA 92008, USA (T) 760.444.4126 (F) 814.238.1608 800-790-2258 (USA & Canada only) www.salimetrics.com support@salimetrics.com</p>
<p>"European Authorized Representative" Biozol Diagnostica Vertrieb GmbH Leipziger Straße 4 85386 Eching, Germany (T) +49 (0)89 3799666-6 www.biozol.de info@biozol.de</p> <p>CE</p>	<p>"United Kingdom Responsible Person" Stratech Scientific Ltd Cambridge House, St Thomas Place, Cambridgeshire Business Park, Ely, CB7 4EX, UK (T) +44 (0) 1638782600 www.stratech.co.uk info@stratech.co.uk</p> <p>UK CA</p>

Updated: August 9th, 2024



101 Innovation Boulevard • Suite 302 • State College, PA 16803
1.800.790.2258 • support@salimetrics.com • salimetrics.com